

# Memtrex\* NY

## pleated filters with hydrophilic Nylon 66 membrane



**figure 1: Memtrex NY filters**

### description and use

Memtrex NY (MNY) filters (figure 1) are uniquely constructed for superior performance with a serial membrane layer design that incorporates a large pore-size membrane layer upstream of the final membrane layer. The first membrane acts like a prefilter to increase the life and efficiency of the cartridge. It is also highly effective in retaining difficult contaminants such as gelatinous particles. The final membrane layer ensures consistent retention characteristics.

MNY filters are designed to ensure maximum downstream cleanliness. SUEZ's Nylon 66 membranes are naturally hydrophilic, non-shedding and do not contain leachable wetting agents. Thermoplastic sealing technologies are used in the cartridge construction eliminating the need for potentially contaminating adhesives. Each cartridge is manufactured under strict production control and is

individually integrity tested. SUEZ is your complete source for filters, crossflow membranes, housings and other filtration equipment.

### typical applications

MNY filters have excellent chemical compatibility making them the ideal choice for a broad range of applications such as high purity water, strong solvents, photoresists and other critical process fluid systems. Typical applications include:

- Filtration of pharmaceutical intermediates
- Filtration of positive photoresists
- Final filtration of beverages

### general properties

Memtrex NY filters are available the following absolute pore size micron ratings: 0.1, 0.2, 0.45, 0.65 and 1.0  $\mu\text{m}$ . Tables 1, 2, 3, 4 and 5 show further details on materials of construction, dimensions, operational limits, integrity testing and flow performance.

**table 1: materials of construction**

description	material of construction
Media	2 layers of Nylon 66 Membrane (prefilter layer over final layer)
Support Layers	Polyester Microfiber
Core and Cage	Polypropylene
Endcaps and Adapters	Polyester

**table 2: dimensions**

nominal O.D.	nominal I.D.	effective filtration area
2.75" (70 mm)	1.25" (31 mm)	6.9 ft <sup>2</sup> (0.64 m <sup>2</sup> )

Find a contact near you by visiting [www.suezwatertechnologies.com](http://www.suezwatertechnologies.com) and clicking on "Contact Us".

\* Trademark of SUEZ; may be registered in one or more countries.

2017 © Copyright SUEZ. All rights reserved.

**table 3: operational limits**


description	operational limits
Maximum Forward Differential Pressure	60 psi (4.1 bar) at 70°F (21°C)
Maximum Reverse Differential Pressure	30 psi (2.1 bar) at 70°F (21°C)
Maximum Operating Temperature	180°F (82°C) at 10 psid (0.7 bar) in water

**table 4: integrity testing**

micron rating	specification
0.1 µm	≤ 50 cc/min at 40 psig (2.76 bar)
0.2 µm	≤ 50 cc/min at 30 psig (2.07 bar)
0.45 µm	≤ 50 cc/min at 15 psig (1.03 bar)
0.65 µm	≤ 50 cc/min at 12 psig (0.83 bar)
1.0 µm	≤ 50 cc/min at 5 psig (0.34 bar)

Air diffusion per 10" module after saturation with clean water

**table 5: flow performance in clean water<sup>1</sup>**

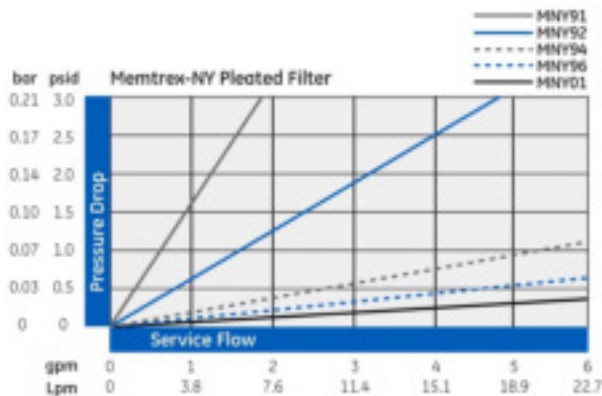

data based on 10" length filter

**table 6: ordering information**

type	micron rating	nominal cartridge length	end #1 adapter	end #2 adapter	elastomer material
MNY	91 = 0.1 µm	1 = 10 in. (25.4 cm)	A = Open End Gasket	A = Open End Gasket	B = Buna-N
	92 = 0.2 µm	2 = 20 in. (50.8 cm)	B = 120 O-Ring	B = 120 O-Ring	E = EPDM
	94 = 0.45 µm	3 = 30 in. (76 cm)	C = 213 O-Ring	C = 213 O-Ring	S = Silicone
	96 = 0.65 µm	4 = 40 in. (101.5 cm)	E = 222 O-Ring	G = Closed End Cap	T = Teflon <sup>2</sup>
	01 = 1.0 µm		F = 226 O-Ring	H = Fin Adapter	Encapsulated (Only in 222 and 226 Sizes)
			J = 020 O-Ring		V = Viton <sup>2</sup>
			Q = 222 O-Ring Stainless Steel Insert		
			Z = 226 O-Ring Stainless Steel Insert		

<sup>1</sup>Teflon and Viton are registered trademarks of DuPont.

**additional information**

- Memtrex NY filters may be autoclaved or in situ steam sterilized (up to 257°F [125°C], 30-minute cycles) for a maximum accumulated exposure of 10 hours. Alternately, the filters may be sanitized with compatible chemical agents.
- SUEZ certifies that the material contained in its Memtrex NY pleated filters meet U.S. FDA requirements for food contact under the applicable regulations in 21 CFR. For further information, contact SUEZ technical services. Memtrex NY filters meet the test criteria for USP class VI-121°C Plastics.
- Aqueous extracts from Memtrex NY filters contain less than 0.25 EU/ml. The filters typically exhibit low levels of non-volatile residues.
- SUEZ filter cartridges are designed and manufactured for resistance to a wide range of chemical solutions. Conditions will vary with each application and users should carefully verify chemical compatibility. Please contact your SUEZ distributor for more information.
- Table 6 provides additional ordering information.

